

## UvA-DARE (Digital Academic Repository)

### Determination of polar 1H-benzotriazoles and benzothiazoles in water by solid-phase extraction and liquid chromatography LTQ FT Orbitrap mass spectrometry

van Leerdam, J.A.; Hogenboom, A.C.; van der Kooi, M.M.E.; de Voogt, P.

**DOI**

[10.1016/j.ijms.2009.02.018](https://doi.org/10.1016/j.ijms.2009.02.018)

**Publication date**

2009

**Published in**

International Journal of Mass Spectrometry

[Link to publication](#)

**Citation for published version (APA):**

van Leerdam, J. A., Hogenboom, A. C., van der Kooi, M. M. E., & de Voogt, P. (2009). Determination of polar 1H-benzotriazoles and benzothiazoles in water by solid-phase extraction and liquid chromatography LTQ FT Orbitrap mass spectrometry. *International Journal of Mass Spectrometry*, 282(3), 99-107. <https://doi.org/10.1016/j.ijms.2009.02.018>

**General rights**

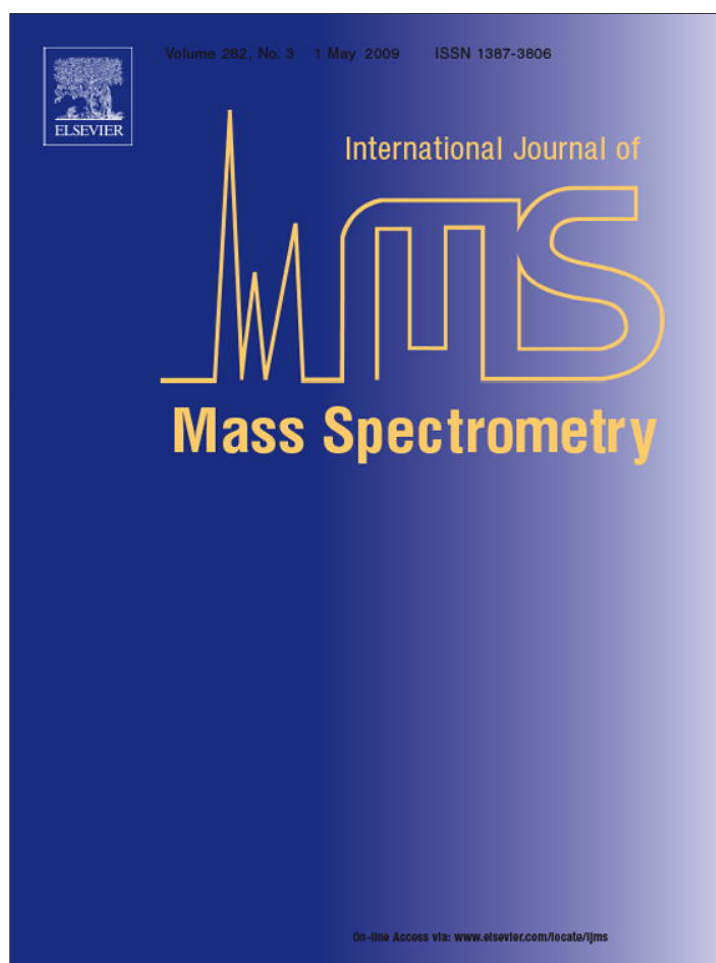
It is not permitted to download or to forward/distribute the text or part of it without the consent of the author(s) and/or copyright holder(s), other than for strictly personal, individual use, unless the work is under an open content license (like Creative Commons).

**Disclaimer/Complaints regulations**

If you believe that digital publication of certain material infringes any of your rights or (privacy) interests, please let the Library know, stating your reasons. In case of a legitimate complaint, the Library will make the material inaccessible and/or remove it from the website. Please Ask the Library: <https://uba.uva.nl/en/contact>, or a letter to: Library of the University of Amsterdam, Secretariat, Singel 425, 1012 WP Amsterdam, The Netherlands. You will be contacted as soon as possible.

*UvA-DARE is a service provided by the library of the University of Amsterdam (<https://dare.uva.nl>)*

Provided for non-commercial research and education use.  
Not for reproduction, distribution or commercial use.



This article appeared in a journal published by Elsevier. The attached copy is furnished to the author for internal non-commercial research and education use, including for instruction at the authors institution and sharing with colleagues.

Other uses, including reproduction and distribution, or selling or licensing copies, or posting to personal, institutional or third party websites are prohibited.

In most cases authors are permitted to post their version of the article (e.g. in Word or Tex form) to their personal website or institutional repository. Authors requiring further information regarding Elsevier's archiving and manuscript policies are encouraged to visit:

<http://www.elsevier.com/copyright>



Contents lists available at ScienceDirect

## International Journal of Mass Spectrometry

journal homepage: [www.elsevier.com/locate/ijms](http://www.elsevier.com/locate/ijms)

# Determination of polar 1H-benzotriazoles and benzothiazoles in water by solid-phase extraction and liquid chromatography LTQ FT Orbitrap mass spectrometry

Jan A. van Leerdam<sup>a,\*</sup>, Ariadne C. Hogenboom<sup>a</sup>, Margaretha M.E. van der Kooi<sup>a</sup>, Pim de Voogt<sup>a,b</sup>

<sup>a</sup> KWR Watercycle Research Institute, P.O. Box 1072, 3430 BB Nieuwegein, The Netherlands

<sup>b</sup> Institute of Biodiversity and Ecosystem Dynamics, University of Amsterdam, Nieuwe Achtergracht 166, 1018 WV Amsterdam, The Netherlands

## ARTICLE INFO

## Article history:

Received 5 November 2008

Received in revised form 10 February 2009

Accepted 10 February 2009

Available online 4 March 2009

## Keywords:

Accurate mass measurement

1H-benzotriazole

Benzothiazole

Linear ion trap (LTQ FT) Orbitrap mass spectrometry

Environmental sample

## ABSTRACT

A sensitive, reliable and robust method for the trace determination of six polar 1H-benzotriazoles and four benzothiazoles in drinking and surface water was developed. These compounds were extracted from water by solid-phase extraction and analyzed by Liquid Chromatography–Electrospray Mass Spectrometry using a linear ion trap–Orbitrap hybrid instrument at high resolution of 30,000 FWHM in the full-scan acquisition mode. At least one product ion was simultaneously detected in the linear ion trap at low mass resolution and was used for confirmation of compound identity.

The compounds studied are soluble in water, resistant to biodegradation, only partially removed in wastewater treatment and they may pass the water treatment processes in the production of drinking water.

The analytes and four internal standards were preconcentrated by solid-phase extraction at low pH. Positive electrospray ionization resulted in protonated molecular ions for all the 1H-benzotriazoles and benzothiazoles.

The mass accuracy was between  $-5$  ppm at  $m/z$  120 and  $-0.1$  ppm at  $m/z$  182 and did not change for more than 2 ppm over a sample sequence of 8 days of analysis time. The optimized method allowed quantifying six benzotriazoles and four benzothiazoles in samples of drinking and surface water down to method detection limits of  $0.01$   $\mu\text{g/L}$ . The recoveries ranged between 45 and 125% in ultrapure, drinking and surface water at a spiking level of  $0.2$   $\mu\text{g/L}$ ; the repeatability was between 2 and 13%. All analytes showed a linear response between  $0.01$  and  $1.0$   $\mu\text{g/L}$ . No significant matrix effect was observed in drinking and surface water, except for the compounds 2-aminobenzothiazole (signal enhancement about 50%) and 2-hydroxybenzothiazole (signal suppression about 25%).

In Dutch drinking water samples, the compounds 1H-benzotriazole, 4- and 5-methyl-1H-benzotriazole, 5,6-dimethyl-1H-benzotriazole, 5-chloro-1H-benzothiazole and benzothiazole were detected. The concentration levels ranged from  $0.01$  to  $0.2$   $\mu\text{g/L}$ . In surface waters, eight out of ten compounds tested were found to be present in concentration levels ranging between  $0.1$  and  $1.0$   $\mu\text{g/L}$ . In addition, in effluents of two sewage treatment plants, eight out of ten compounds tested were present with maximum concentrations for 1H-benzotriazole of  $8$   $\mu\text{g/L}$  and for methyl-1H-benzotriazole of  $3$   $\mu\text{g/L}$  (summed concentration of two isomers).

This work demonstrates the excellent suitability of the LTQ FT Orbitrap mass spectrometer for this type of analysis.

© 2009 Elsevier B.V. All rights reserved.

## 1. Introduction

The benzotriazoles and benzothiazoles are high production volume chemicals that find broad application in various industrial processes as well as in households. The group of 1H-benzotriazoles

are complexing agents that are widely used as anticorrosives, e.g., in engine coolants, aircraft deicers, or anti-freezing liquids, and for silver protection in dish washing liquids. Benzothiazoles are employed in various industrial processes, including rubber manufacturing, and as corrosion inhibitors and fungicides. Accordingly, a release of benzothiazoles into the environment can be expected by industrial discharges from processes where benzothiazoles are involved. Further emissions may stem from the use of benzothiazole-containing consumer products like rubber materials (tires and shoes). The 1H-benzotriazoles and benzothiazoles are soluble in water, resistant

\* Corresponding author.

E-mail address: [ton.van.leerdam@kwrwater.nl](mailto:ton.van.leerdam@kwrwater.nl) (J.A. van Leerdam).

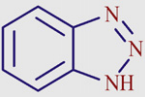
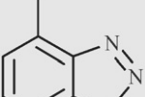
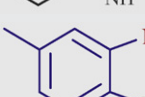
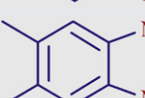
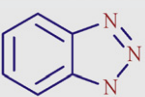
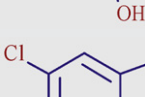
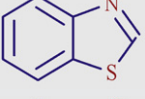
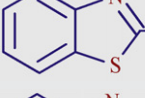
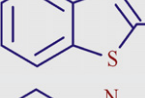
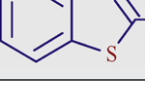
to biodegradation, only partially removed in wastewater treatment [1,2,3,4], and they may pass the water treatment processes employed in the production of drinking water.

Most 1H-benzotriazole and benzothiazole derivatives are polar and thermally labile, which has prevented their detection in the environment as long as gas chromatographic methods have been used. The compounds studied in the present work generally have a  $\log K_{ow}$  below 2.1 (with the exception of methylthio-benzothiazole with a  $\log K_{ow}$  of 3.1, see Table 1). Recently, LC–MS and LC–MS–MS methods have been developed for their measurements in environmental waters. A few reports about the occurrence of 1H-benzotriazoles and benzothiazoles in aquatic compartments have been published since 2003 [6,7]. Giger et al. developed a LC–MS–MS method, using Oasis HLB SPE cartridges for extraction and concentration [5]. Kloepfer et al. [2,3] used a SPE–LC–MS method and showed that benzothiazoles are present in samples from both the influent and effluent of municipal wastewater treatment plants (WWTP) in Berlin and are thus not completely removed

in the wastewater treatment. Weiss and Reemtsma [1] conclude the same for 1H-benzotriazoles found in wastewater and surface water in the Berlin region. The authors expect 1H-benzotriazole and 4-methyl-1H-benzotriazole to be omnipresent contaminants in the water cycle in Germany, due to their limited biodegradability, as high concentrations of several micrograms per liter in treated municipal wastewater were observed.

No information is available until now about the distribution and fate of the 1H-benzotriazoles and benzothiazoles compounds in the Dutch water cycle. Therefore, the objective of this study was to develop a method for the determination of six polar benzotriazoles and four polar benzothiazoles in drinking and surface water on the basis of solid-phase extraction and liquid chromatography–mass spectrometry using a linear ion trap–Orbitrap hybrid instrument. Our main focus was to explore the capabilities of the Orbitrap in the full-scan acquisition mode with high resolution to evaluate if this instrument can achieve sensitivity, accuracy and related performance characteristics required in quantitative anal-

**Table 1**  
Chemical structures, elemental composition, CAS-numbers,  $\log K_{ow}$  and accurate masses of the compounds studied.

Benzotriazoles	Chemical name	Elemental composition	CAS no.	Predicted $\log K_{ow}$ [14]	Accurate mass [M+H] <sup>+</sup>
	1H-benzotriazole	C <sub>6</sub> H <sub>5</sub> N <sub>3</sub>	95-14-7	1.14	120.05562
	4-Methyl-1H-benzotriazole	C <sub>7</sub> H <sub>7</sub> N <sub>3</sub>	29878-31-7	1.60	134.07127
	5-Methyl-1H-benzotriazole	C <sub>7</sub> H <sub>7</sub> N <sub>3</sub>	136-85-6	1.60	134.07127
	5,6-Dimethyl-1H-benzotriazole	C <sub>8</sub> H <sub>9</sub> N <sub>3</sub>	4184-79-6	2.06	148.08692
	1-Hydroxybenzotriazole	C <sub>6</sub> H <sub>5</sub> N <sub>3</sub> O	2592-95-2	0.69	136.05054
	5-Chloro-1H-benzotriazole	C <sub>6</sub> H <sub>4</sub> N <sub>3</sub> Cl	94-97-3	1.74	154.01665
<b>Benzothiazoles</b>					
	Benzothiazole	C <sub>7</sub> H <sub>5</sub> NS	95-16-9	2.01	136.02155
	2-Aminobenzothiazole	C <sub>7</sub> H <sub>6</sub> N <sub>2</sub> S	136-95-8	1.89	151.03245
	2-Hydroxybenzothiazole	C <sub>7</sub> H <sub>5</sub> NSO	934-34-9	1.76	152.01646
	2-(Methylthio)benzo-thiazole	C <sub>8</sub> H <sub>7</sub> NS <sub>2</sub>	615-22-5	3.10	182.00927

ysis of low level water contaminants. Nominal mass product ions were used for confirmation of the identity of the compounds.

This paper discusses the performance of the quantitative method for the determination of six 1H-benzotriazoles and four benzothiazoles in surface water and drinking water. The benzotriazoles selected include 1H-benzotriazole, 4-methyl-1H-benzotriazole, 5-methyl-1H-benzotriazole, 5,6-dimethylbenzotriazole and 5-chloro-1H-benzotriazole. Benzothiazole, 2-hydroxybenzothiazole, 2-aminobenzothiazole and 2-(methylthio)benzothiazole are the benzothiazoles that were investigated. Table 1 presents the structures of the compounds. The compound 2-mercaptobenzothiazole was not studied due to its oxidation in aqueous solutions [3].

The performance of the overall off-line SPE and LC–MS method was evaluated in terms of linearity, limit of detection (LOD) and quantification (LOQ), repeatability and recovery. In addition, to obtain an optimum recovery for all benzotriazoles and benzothiazoles, different SPE materials (type, amount and volume), elution solvents, evaporation temperatures for concentration the sample extract and the influence of the pH of the water sample were evaluated. Hereto, water samples, i.e., ultrapure water, drinking water and surface water were spiked with the analytes at different concentration levels. The influence of the matrix effect in the quantitative analytical LC–Orbitrap MS method was also studied.

Finally, the optimized method was used to monitor the presence of the benzotriazoles and benzothiazoles in samples of drinking water, surface waters and WWTP effluents.

## 2. Experimental

### 2.1. Reagents and chemicals

1H-benzotriazole was purchased from Riedel-de-Haën (Sigma–Aldrich, Zwijndrecht, The Netherlands) and Acros (Boom, Meppel, The Netherlands) supplied 5-methylbenzotriazole, 5,6-dimethylbenzotriazole and benzothiazole. 5-chloro-1H-benzotriazole and 2-hydroxybenzothiazole were obtained from Aldrich (Sigma–Aldrich, Zwijndrecht, The Netherlands).

2-Aminobenzothiazole and 2-(methylthio)benzothiazole were supplied by Fluka (Sigma–Aldrich, Zwijndrecht, The Netherlands).

Two compounds, 4-methylbenzotriazole and the internal standard (I.S.) 1H-benzotriazole- $d_4$ , are not commercially available. Therefore they were synthesized by Mercachem (Nijmegen, The Netherlands). The reference standards were purchased in the highest purity available (96% or greater). The I.S. and a second one, the compound atrazine- $d_5$  (CDN, J.H. Ritmeester B.V., Utrecht, The Netherlands) were added to the sample extract to check the performance of the LC–MS analysis. 1H-benzotriazole- $d_4$ , was used for correction of the calculated concentrations.

The compounds fenuron and chloroxuron (Dr. Ehrenstorfer, Boom, Meppel, The Netherlands) were used as internal standards to monitor the sample preparation process and were added to the water samples. Water used for LC–MS analysis was generated from an ultrapure water system from Millipore (Bedford, MA, USA) with a specific resistance of 18 M $\Omega$  cm. Ultra gradient acetonitrile (Mallinckrodt Baker, Deventer, The Netherlands) was used as an organic modifier.

Stock solutions of the individual compounds were prepared in acetonitrile (200 mg/L) and working solutions (0.01–1.0 mg/L) were prepared in a mixture of ultrapure water and acetonitrile (50:50; v/v). Stock solutions were stored at a temperature of  $-18^\circ\text{C}$  and working solutions at  $4^\circ\text{C}$ . The ultrapure water was obtained by purifying demineralized water in a Milli-Q-system. All organic solvents were HPLC-grade. A 30% hydrochloric acid-solution was used

and the quality of the formic acid was 98%. The sea sand had been glowd and ignited by the manufacturer.

### 2.2. Solid-phase extraction

Preconcentration of the water samples was performed by comparing two different commercially available cartridges with adsorbent: Oasis HLB (5 mL glass cartridges, 200 mg) and Oasis HLB (6 mL polypropylene, 500 mg; Waters, Milford, MA, USA). A 5 mL glass or 6 mL polypropylene cartridge, filled with 4 gram of sea sand (Mallinckrodt Baker, Deventer, The Netherlands), was placed above the cartridge with adsorbent, in order to prevent the SPE material from clogging. The SPE cartridges were conditioned successively with 10 mL of acetonitrile, 10 mL of methanol and 10 mL of ultrapure water (pH 2, except for the neutral pH experiments). The sea sand cartridges were conditioned in the same way, but the methanol was omitted, due to the absence of active sites. In order to examine the influence of the pH on the recoveries, spiked samples of drinking water and surface water were concentrated at a pH between 7.8 and 8.0, and at pH 2.

The internal standard compounds fenuron and chloroxuron were spiked to 1 L of the water sample at a concentration level of 1  $\mu\text{g/L}$ . The samples were acidified to pH 2 with 30% hydrochloric acid (except for the experiments at neutral pH) and the water samples were percolated through the cartridges with a flow of 10 mL/min. After sampling the water, the sea sand cartridges were removed and the SPE cartridges were washed twice with 5 mL of acidified water (pH 2). The cartridges were dried under vacuum for 90 min until dryness and eluted 3 times with 2.5 mL of organic solvent. Three different solvents were evaluated: acetonitrile, methanol and a mixture of acetonitrile/methanol (50:50; v/v). The eluates were concentrated to a final volume of 500  $\mu\text{L}$  under a gentle flow of nitrogen. The temperature during the concentration step was evaluated (37, 45 and  $56^\circ\text{C}$ ). Lastly, the extracts were diluted with 500  $\mu\text{L}$  internal standard solution (benzotriazole- $d_4$  and atrazine- $d_5$ ) in ultrapure water at a concentration level of 2 mg/L. The final concentration of the deuterated internal standards in the sample extract is 1 mg/L.

### 2.3. Chromatographic separation

The liquid chromatograph was equipped with a Surveyor autosampler model Plus and a Surveyor quaternary gradient HPLC-pump (Thermo Fisher Scientific, Breda, The Netherlands). The separation of the compounds was performed using a 250 mm  $\times$  4.6 mm i.d. Phenomenex Aqua  $C_{18}$ -column with 5  $\mu\text{m}$  particles. The precolumn used was a 4.0  $\times$  3.0 mm i.d. Phenomenex Security Guard column (both columns were supplied by Bester, Amsterdam, The Netherlands). These columns were conditioned in a column thermostat that was maintained at a temperature of  $21^\circ\text{C}$ . 20  $\mu\text{L}$  of the water/acetonitrile-extract was injected into the LC system and the analytes were separated using a linear gradient elution at a flow rate of 1.0 mL/min. The gradient started with 5% acetonitrile/95% water/0.05% formic acid, increased to 40% acetonitrile/60% water/0.05% formic acid in 15 min, then increased to 100% acetonitrile in 10 min, and was held constant for 10 min. Between the runs, the analytical column was re-equilibrated for 15 min.

### 2.4. Mass spectrometry

A linear ion trap-Fourier Transform (LTQ-FT) Orbitrap mass spectrometer (Thermo Electron, Bremen, Germany) was used. The ion trap part of this system was equipped with an Ion Max Elec-



rospray ionization (ESI) probe that was applied in the positive-ion mode for these compounds.

From the LC column effluent, 200  $\mu\text{L}/\text{min}$  was introduced into the mass spectrometer, whereas the remaining 800  $\mu\text{L}/\text{min}$  was directed to waste. The full-scan accurate mass spectra from 100 to 1000 Da, which were obtained at a resolution of 30,000 FWHM, were processed using the Xcalibur version 2.0 software. The total cycle time depends upon the resolution; at a resolution of 30,000 FWHM the total cycle time is about 0.55 s.

Before measuring the samples, accurate mass calibration was performed using flow injection of a 1,3,6-polytyrosine-solution ( $m/z$  181, 507 and 997) in methanol/water (50/50; v/v) with 0.1% formic acid at a flow rate of 10  $\mu\text{L}/\text{min}$ .

The conditions in the Electrospray positive-ion mode were: source voltage 3.6 kV, heated capillary temperature 275  $^{\circ}\text{C}$ , capillary voltage 30 V and tube lens 70 V.

In the LTQ component of the instrument, the FT part was set to 26  $^{\circ}\text{C}$  and helium was used as damping gas. All measurements were done using the automatic gain control (AGC) of the LTQ to adjust the number of ions entering the trap. Product ions were generated in the LTQ trap at a normalized collision energy setting of 35% and using an isolation width of 2 Da.

In order to investigate the effect of the matrix on the response of the mass spectrometer, the signal suppression was calculated for all compounds in sample extracts of ultrapure water, drinking water and surface water compared to a standard in ultrapure water.

Identification of the compounds was performed with the accurate mass of the protonated molecule within a mass window of 7 ppm together with one product ion (nominal mass). The retention times of the compounds were compared to those of the compounds in the calibration standard solution of the final analysis.

Quantification of the compounds was obtained by using a 7-points external calibration curve ranging between 0.01 and 1.0 mg/L, corresponding to 0.01–1.0  $\mu\text{g}/\text{L}$  water. The solutions for the calibration curve were prepared in acetonitrile/water (50/50; v/v). Peak areas of the external calibration curve were corrected for the internal standard benzotriazole- $d_4$ . The recoveries of the internal standards fenuron and chloroxuron were used to check the performance of the extraction procedure, whereas atrazine- $d_5$  was measured to check the performance of the LC-MS analysis. Table 1 shows the chemical structure, chemical name, elemental composition, CAS number, predicted  $\log K_{ow}$  value [14] and the accurate masses of the protonated compounds.

## 2.5. Method evaluation

The performance of the method was evaluated in terms of linearity, recovery, limit of detection, limit of quantification and

repeatability. In addition, matrix effects were studied and quantified. Furthermore, confirmation of the target compounds was assessed based on concept of the use of identification point criteria [10].

The optimized method was used for the analyses of different drinking waters, surface waters and WWTP effluents originating from different locations in The Netherlands.

## 3. Results and discussion

### 3.1. Extraction optimization

In order to examine the influence of the pH on the recovery, spiked samples of drinking water and surface waters were concentrated at neutral pH (between 7.8 and 8.0) and pH 2. Besides the two pH values, two spiking levels (0.1 and 1.0  $\mu\text{g}/\text{L}$ ) were examined. The samples were processed, as described earlier in Section 2, using 6 mL glass cartridges packed with 200 mg Oasis HLB and acetonitrile as elution solvent. The recovery of all benzotriazoles at neutral pH was less than 50%, whereas the recovery at pH 2 was in the range 50–110%. For the benzothiazoles studied, the recovery of only 2-aminobenzothiazole decreases from >50% to <30% when lowering the pH to 2, due to the basic properties of this compound. The recovery of the other benzothiazoles was not significantly affected by the pH (difference less than 5%).

The average recovery of all the compounds (both benzotriazoles and benzothiazoles) extracted at neutral pH is 42%, whereas the average recovery at pH 2 is 65%. In further experiments, water samples were acidified to pH 2.

In the second set of experiments, the effects of the amount of adsorbent and the solvent composition for the elution of the SPE-cartridges were investigated. Therefore, two different cartridge volumes of the Oasis HLB type and three different compositions of elution solvents were evaluated. Drinking water, spiked with 0.1  $\mu\text{g}/\text{L}$  of the target compounds, was processed as described in Section 2. Besides the use of cartridges that were filled with 200 and 500 mg Oasis HLB, the elution solvents acetonitrile, a mixture of acetonitrile/methanol (50:50; v/v) and methanol were assessed. All experiments were performed in duplicate and the average recoveries are shown in Table 2. The results shows that target compound recoveries are higher than 68% when using 500 mg adsorbent and the acetonitrile/methanol elution solvent. Every other combination causes a recovery below 34% for one or more compounds. Therefore, 500 mg adsorbent and a mixture of acetonitrile/methanol (50:50; v/v) were used in the further evaluation of the performance of the method.

To study the possible loss of compounds during the evaporation step of the sample extracts, three different evaporation temperatures were investigated (37, 45 and 56  $^{\circ}\text{C}$ ). To that end, 7.5 mL of a

**Table 2**  
Effect of adsorbent amount and composition of the elution solvent on the mean analyte recovery. Spike level 0.1  $\mu\text{g}/\text{L}$  in drinking water ( $n = 2$ ).

Compound name	200 mg/methanol	200 mg/acetonitrile: methanol 50:50	200 mg/acetonitrile	500 mg/methanol	500 mg/acetonitrile: methanol 50:50	500 mg/acetonitrile
1H-benzotriazole	74 $\pm$ 4	74 $\pm$ 1	68 $\pm$ 0	73 $\pm$ 0	68 $\pm$ 3	41 $\pm$ 1
4-Methyl-1H-benzotriazole	87 $\pm$ 0	88 $\pm$ 4	79 $\pm$ 1	81 $\pm$ 1	79 $\pm$ 4	55 $\pm$ 6
5-Methyl-1H-benzotriazole	78 $\pm$ 1	76 $\pm$ 4	68 $\pm$ 1	71 $\pm$ 2	70 $\pm$ 0	38 $\pm$ 1
5,6-Dimethyl-1H-benzotriazole	90 $\pm$ 0	88 $\pm$ 2	64 $\pm$ 13	85 $\pm$ 1	82 $\pm$ 3	28 $\pm$ 3
5-Chloro-1H-benzotriazole	77 $\pm$ 2	77 $\pm$ 2	59 $\pm$ 12	73 $\pm$ 0	72 $\pm$ 1	13 $\pm$ 1
Benzothiazole	91 $\pm$ 8	86 $\pm$ 4	87 $\pm$ 11	116 $\pm$ 3	132 $\pm$ 18	116 $\pm$ 10
2-Aminobenzothiazole	32 $\pm$ 0	32 $\pm$ 2	21 $\pm$ 2	69 $\pm$ 0	73 $\pm$ 5	13 $\pm$ 7
2-Hydroxybenzothiazole	78 $\pm$ 4	76 $\pm$ 4	79 $\pm$ 15	75 $\pm$ 5	85 $\pm$ 0	61 $\pm$ 5
2-(Methylthio)benzothiazole	89 $\pm$ 1	95 $\pm$ 3	70 $\pm$ 11	34 $\pm$ 0	98 $\pm$ 5	81 $\pm$ 3
Fenuron	87 $\pm$ 2	84 $\pm$ 2	99 $\pm$ 12	86 $\pm$ 3	95 $\pm$ 4	104 $\pm$ 9
Chloroxuron	119 $\pm$ 4	118 $\pm$ 2	104 $\pm$ 17	113 $\pm$ 3	120 $\pm$ 3	109 $\pm$ 0
Atrazine- $d_5$	114 $\pm$ 2	110 $\pm$ 2	95 $\pm$ 18	110 $\pm$ 2	116 $\pm$ 4	94 $\pm$ 6

**Table 3**

Retention time and product ion masses measured in LTQ with intensity at a normalized collision energy of 35, 40 and 45%.

Chemical name	Rt (min)	LTQ mass 1	LTQ mass 2	Ion ratio mass 2/mass 1 CE = 35%	Ion ratio mass 2/mass 1 CE = 40%	Ion ratio mass 2/mass 1 CE = 45%
1H-benzotriazole ( <i>m/z</i> 120)	15.0	120	92	5	15	35
4-Methyl-1H-benzotriazole ( <i>m/z</i> 134)	18.9	134	106	5	10	15
5-Methyl-1H-benzotriazole ( <i>m/z</i> 134)	19.1	134	106	5	10	15
5,6-Dimethyl-1H-benzotriazole ( <i>m/z</i> 148)	22.1	148	120	15	40	40
1-Hydroxy-benzotriazole ( <i>m/z</i> 136)	11.0	136	119	5	20	50
5-Chloro-1H-benzotriazole ( <i>m/z</i> 154)	21.9	154	126	15	50	80
Benzothiazole ( <i>m/z</i> 136)	24.7	136	109	0	10	50
2-Aminobenzothiazole ( <i>m/z</i> 151)	11.0	151	124	10	60	60
2-Hydroxybenzothiazole ( <i>m/z</i> 152)	22.7	152	124	5	10	15
2-(Methylthio)benzotriazole ( <i>m/z</i> 182)	28.4	182	167	20	100	100

standard solution of all the compounds in acetonitrile at a concentration level of 0.1 mg/L was evaporated to a final volume of about 0.5 mL. The recovery of all compounds was in the range of 95–100% related to a standard solution of 0.1 mg/L in acetonitrile. In conclusions the temperature of the evaporation step did not influence the recovery of the analytes. For practical reasons, a temperature of 45 °C for the evaporation step was selected.

### 3.2. Method performance

A chromatographic separation of the isomers 4-methyl-1H-benzotriazole and 5-methyl-1H-benzotriazole was achieved with the chromatographic conditions used (see Section 2). The individual detection of these two isomers is important due to their differences in toxicity and degradation behavior [13].

#### 3.2.1. Confirmation

For confirmation of target compounds, LC relative retention time criteria (retention time window <2.5%) and mass spectrometric identification criteria need to be fulfilled. The latter being based on the concept of identification points [10]. For accurate mass screening using ToF or Orbitrap MS instruments, no criteria are described and recently some proposals were made by Nielen et al. [11]. For high resolution screening (resolution  $\geq 20,000$  and a mass accuracy  $\leq 5$  mDa) they proposed 2 identification points. Each product ion (low mass resolution MS) contributes also 2 points. Thus acquiring a high resolution precursor ion in combination with at least one product ion and the LC relative retention time meets the minimum requirement of 4 identification points. In Table 3, the two most intense (nominal mass) ions generated in the LTQ part of the instrument by normalized collision energy of 35, 40 and 45% were tabulated. At a default normalized collision energy value of 35%, the precursor ion is the most dominant ion beside some product ions. One product ion was used to check the identity of the compounds. Table 4 presents the accurate mass of the measured product ion with the elemental composition. For all benzotriazoles

except 1-hydroxybenzotriazole, a loss of N<sub>2</sub> was observed. For the benzothiazoles, losses of CHN, CO and SCH<sub>3</sub> were observed.

### 3.3. Matrix effects

The matrix effects (the possibility of signal suppression or enhancement) were assessed using a method described by Matuszewski et al. [9].

To that end, the peak area of the analyte in the sample extract spiked after extraction (B) was compared with the peak area of the corresponding analyte in the standard solution (A). The matrix effects (ME) were calculated as follows:

$$ME(\%) = \frac{B}{A \times 100}$$

Consequently, a value higher than 100% means that signal enhancement occurs in the mass spectrometer due to the co-extracted and co-eluted matrix. A value lower than 100% means that signal suppression occurs in the mass spectrometer. This matrix effect was measured in twofold and calculated for ultrapure water, drinking water and surface water at a concentration level of 0.2 mg/L for each compound. The results were corrected for the blank values. Table 5 shows the results of these experiments.

No significant matrix effect was observed in drinking and surface water, except for the compounds 2-aminobenzothiazole (signal enhancement about 50%) and 2-hydroxybenzothiazole (signal suppression about 25%).

#### 3.3.1. Overall recovery and precision

The recovery and precision of the overall method were evaluated by analyzing extracts of ultrapure, drinking water and surface water by LC-Orbitrap MS, spiked before extraction with 0.025 and 0.2 and 1.0  $\mu\text{g/L}$ , respectively, of the analytes. For the 0.025 and 0.2  $\mu\text{g/L}$  level, seven replicates were analyzed. For the 1.0  $\mu\text{g/L}$  spiking level, only two replicates were analyzed, therefore the precision at that level is less accurate.

**Table 4**

Accurate masses and elemental composition of product ion masses measured.

Chemical name	Elemental composition of [M+H] <sup>+</sup>	Accurate mass of product ion	Elemental composition of product ion	Mass error (ppm)	Loss of
1H-benzotriazole ( <i>m/z</i> 120)	C <sub>6</sub> H <sub>6</sub> N <sub>3</sub>	92.04874	C <sub>6</sub> H <sub>6</sub> N	-8.0	N <sub>2</sub>
4-Methyl-1H-benzotriazole ( <i>m/z</i> 134)	C <sub>7</sub> H <sub>8</sub> N <sub>3</sub>	106.06471	C <sub>7</sub> H <sub>8</sub> N	-3.9	N <sub>2</sub>
5-Methyl-1H-benzotriazole ( <i>m/z</i> 134)	C <sub>7</sub> H <sub>8</sub> N <sub>3</sub>	106.06471	C <sub>7</sub> H <sub>8</sub> N	-3.9	N <sub>2</sub>
5,6-Dimethyl-1H-benzotriazole ( <i>m/z</i> 148)	C <sub>8</sub> H <sub>10</sub> N <sub>3</sub>	120.08041	C <sub>8</sub> H <sub>10</sub> N	-3.0	N <sub>2</sub>
1-Hydroxy-benzotriazole ( <i>m/z</i> 136)	C <sub>6</sub> H <sub>6</sub> N <sub>3</sub> O	119.04742	C <sub>6</sub> H <sub>5</sub> N <sub>3</sub>	-3.2	OH
5-Chloro-1H-benzotriazole ( <i>m/z</i> 154)	C <sub>6</sub> H <sub>5</sub> N <sub>3</sub> Cl	126.01014	C <sub>6</sub> H <sub>5</sub> NCl	-2.9	N <sub>2</sub>
Benzothiazole ( <i>m/z</i> 136)	C <sub>7</sub> H <sub>6</sub> NS	109.01012	C <sub>7</sub> H <sub>5</sub> S	-4.8	CHN
2-Aminobenzothiazole ( <i>m/z</i> 151)	C <sub>7</sub> H <sub>7</sub> N <sub>2</sub> S	124.02114	C <sub>6</sub> H <sub>6</sub> NS	-3.2	CHN
2-Hydroxybenzothiazole ( <i>m/z</i> 152)	C <sub>7</sub> H <sub>6</sub> NSO	124.02114	C <sub>6</sub> H <sub>6</sub> NS	-3.2	CO
2-(Methylthio)benzotriazole ( <i>m/z</i> 182)	C <sub>8</sub> H <sub>8</sub> NS <sub>2</sub>	166.98579	C <sub>7</sub> H <sub>5</sub> NS	0.0	SCH <sub>3</sub>

**Table 5**

Matrix effect (ME; %) according to Matuszewski et al.[9] for target compounds and internal standards in different types of water. Concentration tested 1.0 µg/L; *n* = 2. ME < 100%: signal suppression, ME > 100%: signal enhancement.

Compound name	Ultrapure water	Drinking water	Surface water
1H-benzotriazole	103	118	95 <sup>a</sup>
4-Methyl-1H-benzotriazole	102	126	110 <sup>a</sup>
5-Methyl-1H-benzotriazole	95	109	108 <sup>a</sup>
5,6-Dimethyl-1H-benzotriazole	108	118	114 <sup>a</sup>
1-Hydroxybenzotriazole	87	103	106
5-Chloro-1H-benzotriazole	89	99	101 <sup>a</sup>
Benzothiazole	100	114	106 <sup>a</sup>
2-Aminobenzothiazole	132	153	143
2-Hydroxybenzothiazole	91	74	78 <sup>a</sup>
2-(Methylthio)benzothiazole	100	98	101 <sup>a</sup>
Fenuron	76	70	66
Chloroxuron	101	100	96
1H-benzotriazole-d <sub>4</sub>	85	103	106
Atrazine-d <sub>5</sub>	88	83	84

<sup>a</sup> Corrected for blank value.

After isolation with Oasis HLB, the extracts were measured in duplicate and the recoveries were calculated and corrected for the I.S. 1H-benzotriazole-d<sub>4</sub>. Table 6 shows the average recoveries at three concentration levels together with the standard deviation (RSD). No correction for the matrix effect is performed.

The overall recovery numbers are equal to or below the matrix effect numbers, indicating (1) that the matrix effect is an important contribution to the overall recoveries and (2) that losses of analytes do also occur during sample extraction and treatment.

At a concentration level of 0.2 µg/L, the recovery in ultrapure water ranged from 45% for 2-(methylthio)-benzothiazole up to 97% for the compound 5,6-dimethyl-1H-benzotriazole. In drinking water, the recoveries were in general slightly higher than in ultrapure water. In surface water, the recoveries ranged from 57% for 1-hydroxybenzotriazole up to 125% for 4-methyl-1H-benzotriazole. The relatively high recovery of 4-methyl-benzotriazole (125%) can be explained by signal enhancement (ME = 110% in surface water at a spiking level of 1.0 µg/L) due to matrix effects.

The reproducibility of the total procedure was good with relative standard deviations of 3–7% in ultrapure water, 4–8% in drinking water and 2–13% in surface water at a concentration level of 0.2 µg/L.

### 3.4. Linearity

A 7-point standard calibration curve in water/acetonitrile (50/50) was used for the quantification of the compounds. In order to evaluate the linear range, the *R*<sup>2</sup> of the concentration–response

**Table 6**

Recoveries (%) and RSD (between brackets, %) of the overall analytical method (*n* = 7).

Compound name	Ultrapure water			Drinking water			Surface water		
	0.025 µg/L	0.2 µg/L	1.0 µg/L <sup>a</sup>	0.025 µg/L	0.2 µg/L	1.0 µg/L <sup>a</sup>	0.025 µg/L	0.2 µg/L	1.0 µg/L <sup>a</sup>
1H-benzotriazole	38(11)	75(4)	81	54(12)	86(5)	91	256 <sup>b</sup> (47)	115 <sup>b</sup> (13)	88 <sup>b</sup>
4-Methyl-1H-benzotriazole	74(5)	76(4)	90	73(3)	90(3)	116	193 <sup>b</sup> (14)	125 <sup>b</sup> (4)	104 <sup>b</sup>
5-Methyl-1H-benzotriazole	54(7)	81(3)	89	38(13)	87(4)	94	94 <sup>b</sup> (32)	99 <sup>b</sup> (6)	91 <sup>b</sup>
5,6-Dimethyl-1H-benzotriazole	86(4)	97(2)	100	57(5)	91(4)	94	57(5)	88(2)	87 <sup>b</sup>
1-Hydroxybenzotriazole	99(3)	68(2)	67	49(5)	55(3)	63	88(5)	57(3)	67
5-Chloro-1H-benzotriazole	84(4)	84(3)	107	68(3)	75(3)	108	72 <sup>b</sup> (4)	76 <sup>b</sup> (5)	98 <sup>b</sup>
Benzothiazole	71(6)	62(3)	65	102(10)	93(6)	87	101 <sup>b</sup> (13)	111 <sup>b</sup> (2)	107 <sup>b</sup>
2-Aminobenzothiazole	61(4)	77(3)	61	50(6)	68(4)	62	75(7)	75(3)	72
2-Hydroxybenzothiazole	70(4)	76(2)	89	54(7)	62(4)	74	79 <sup>b</sup> (9)	67 <sup>b</sup> (4)	69 <sup>b</sup>
2-(Methylthio)benzothiazole	24(5)	45(7)	52	58(7)	74(8)	82 <sup>b</sup>	76 <sup>b</sup> (9)	94 <sup>b</sup> (10)	67 <sup>b</sup>

<sup>a</sup> *n* = 2.

<sup>b</sup> Corrected for the blank value.

curve was calculated for the analytes in ultrapure water and in extracts of ultrapure, drinking water and surface water, spiked before extraction. In total, 10 analytes in 4 matrices were analyzed. The concentrations ranged between 0.01 and 1.0 µg/L (sample). The *R*<sup>2</sup> ranges invariably between 0.992 and 1.000. The lowest value for the *R*<sup>2</sup> was observed for the compound 5-chloro-1H-benzotriazole (averaged value over the 4 matrices is 0.994). In the majority of cases (31 of 40), the *R*<sup>2</sup> was higher than 0.996. This means that the optimized LC–MS method can be used for the quantitative analyses of the analytes in the range 0.01–1.0 µg/L in both drinking water and surface water.

### 3.4.1. LOD and LOQ

Quantification and detection limits were determined using a signal-to-noise approach. The limit of quantification was estimated for the concentration of compound that results in a signal-to-noise ratio of 10:1. The limit of detection corresponds to the concentration that results in a signal-to-noise ratio of 3:1.

Hereto, a standard dissolved in ultrapure water, and extracts of ultrapure water, drinking water and surface water, spiked after extraction with 0.025 mg/L of the analytes, were analyzed 7 times. The signal-to-noise ratios were measured and the average values are presented in Table 7. The data in this table show that the limit of detection ranges between <1 and 3 ng/L for ultrapure water, from <1 to 8 ng/L for drinking water and from <1 to 10 ng/L in surface waters. The limit of quantification ranges between 1 and 17 ng/L for ultrapure water, 1 and 9 ng/L for drinking water and between 1 and 33 ng/L for surface waters.

### 3.5. Mass accuracy

While for identification of unknowns the mass accuracy is of largest importance [8], this is not a strict prerequisite for target analysis such as in this study, as the mass error could be accounted for from the injection of calibration solutions. Still, the mass accuracy should be stable over a whole analytical sequence to ensure accurate peak detection and integration. To study the mass accuracy, the drift in mass error of the target compounds with low molecular mass over a time period of 8 days was evaluated. The accurate mass was derived from the apex of the chromatographic peak. The accurate mass at the start of the experiments and after 3 and 8 days are presented in Table 8. The concentration of the compounds is 0.1 mg/L. The mass spectrometer was calibrated externally (see Section 2) before the start of the experiments. No recalibration was performed during the experiments.

In general, the mass error for all the compounds ranged from –0.1 to –4.8 ppm for all studied ions and was largest at the lowest *m/z* range. The mass accuracy of the compound with the lowest molecular mass (1H-benzotriazole, *m/z* 120) was –4.8 ppm. After



**Table 7**  
Limits of detection (LOD) and limits of quantitation (LOQ) in standards and extracts.

Compound	Standard		Extract of spiked ultrapure water	
	LOD <sup>a</sup> (ng/L)	LOQ <sup>b</sup> (ng/L)	LOD <sup>a</sup> (ng/L)	LOQ <sup>b</sup> (ng/L)
1H-benzotriazole	3	10	3	10
4-Methyl-1H-benzotriazole	<1	2	<1	2
5-Methyl-1H-benzotriazole	<1	2	<1	2
5,6-Methyl-1H-benzotriazole	<1	1	<1	1
1-Hydroxybenzotriazole	1	4	<1	3
5-Chloro-1H-benzotriazole	<1	2	<1	2
Benzothiazole	3	10	5	17
2-Aminobenzothiazole	4	13	1	4
2-Hydroxy-benzothiazole	3	10	3	9
2-(Methylthio)benzothiazole	<1	3	<1	3

Compound	Drinking water		Surface water	
	LOD <sup>a</sup> (ng/L)	LOQ <sup>b</sup> (ng/L)	LOD <sup>a</sup> (ng/L)	LOQ <sup>b</sup> (ng/L)
1H-benzotriazole	4	9	5	11
4-Methyl-1H-benzotriazole	1	4	1	2
5-Methyl-1H-benzotriazole	1	4	1	4
5,6-Methyl-1H-benzotriazole	<1	2	2	7
1-Hydroxybenzotriazole	<1	3	2	5
5-Chloro-1H-benzotriazole	2	6	2	7
Benzothiazole	6	21	10	33
2-Aminobenzothiazole	1	5	2	5
2-Hydroxy-benzothiazole	8	25	6	20
2-(Methylthio)benzothiazole	<1	1	<1	1

<sup>a</sup> s/n = 3.<sup>b</sup> s/n = 10.

8 days there was a drift in mass accuracy of  $-1.7$  ppm. The mass accuracy of the compound with the highest molecular mass (2-(methylthio)benzothiazole ( $m/z$  182)) was  $-0.1$  ppm. Typically, the mass error did not change more than  $-2$  ppm over the course of a sequence of more than 250 injections corresponding to a total analysis time of 8 days. Any effect of the signal intensity on mass accuracy was not noticed, indicating the large dynamic range of the Orbitrap MS with respect to mass accuracy. Similar results were obtained by Krauss and Hollender [12], who used the Orbitrap for the trace level determination of low molecular weight nitrosamines.

On the basis of these results, peak integration for quantitative analysis was done from ion chromatograms extracted for each ion at a range of  $\pm 7$  ppm around the theoretical  $m/z$  value. If a better mass accuracy at the low  $m/z$  range would be needed, a lower calibration mass would be required than that of  $m/z$  181 in the routine calibration solution.

### 3.6. Application to real water samples

Ten samples of drinking water were collected in 2007 from ten different locations. Also ten different surface water

locations and effluents from two sewage treatment plants, originating from all over The Netherlands, were sampled and analyzed with the method developed. All samples were stored at  $4^{\circ}\text{C}$  for a maximum period of 7 days in green glass bottles prior to analysis. In drinking water samples (originating from surface waters) from The Netherlands, the compounds 1H-benzotriazole, 4- and 5-methyl-1H-benzotriazole, 5,6-dimethyl-1H-benzotriazole, 5-chloro-1H-benzotriazole and benzothiazole were detected (see Fig. 1). The concentration levels ranged from 0.01 for the compounds 5-methyl-1H-benzotriazole, 5,6-dimethyl-1H-benzotriazole, 5-chloro-1H-benzotriazole and benzothiazole to  $0.2 \mu\text{g/L}$  for 1H-benzotriazole and 4-methyl-1H-benzotriazole. In particular, the latter two compounds have the potential to reach drinking water prepared either directly from surface water or from surface water via bank filtration and should, therefore, be included in future monitoring studies. These results are in accordance with a study on occurrence of benzotriazoles and benzothiazoles in different samples from the Berlin region, reported by Weiss and Reemtsma [1], who also found the highest concentrations for 1H-benzotriazole and 4-methyl-1H-benzotriazole. The concentration of the isomer 5-methyl-1H-benzotriazole is significantly lower due to its lower microbial stability [1].

**Table 8**  
Drift in mass error (in ppm) within an analysis time range of 8 days without mass calibration. Compounds are presented in decreasing order of mass error. Concentration 0.1 mg/L.

Compound name	Mass error (ppm), start	Mass error (ppm), 3 days	Mass error (ppm), 8 days	Total drift (ppm)
1H-benzotriazole ( $m/z$ 120)	-4.8	-5.3	-6.5	-1.7
4-Methyl-1H-benzotriazole ( $m/z$ 134)	-4.0	-3.3	-4.8	-0.9
5-Methyl-1H-benzotriazole ( $m/z$ 134)	-2.8	-4.4	-4.8	-2.0
1-Hydroxybenzotriazole ( $m/z$ 136)	-2.8	-2.5	-3.7	-1.0
Benzothiazole ( $m/z$ 136)	-2.6	-2.3	-4.0	-1.4
5-Chloro-1H-benzotriazole ( $m/z$ 154)	-2.1	-2.4	-3.6	-1.5
2-Aminobenzothiazole ( $m/z$ 151)	-2.1	-2.5	-3.3	-1.3
2-Hydroxybenzothiazole ( $m/z$ 152)	-1.9	-2.0	-3.0	-1.1
5,6-Dimethyl-1H-benzotriazole ( $m/z$ 148)	-1.9	-2.4	-3.6	-1.8
2-(Methylthio)benzothiazole ( $m/z$ 182)	-0.1	-0.4	-1.8	-1.7

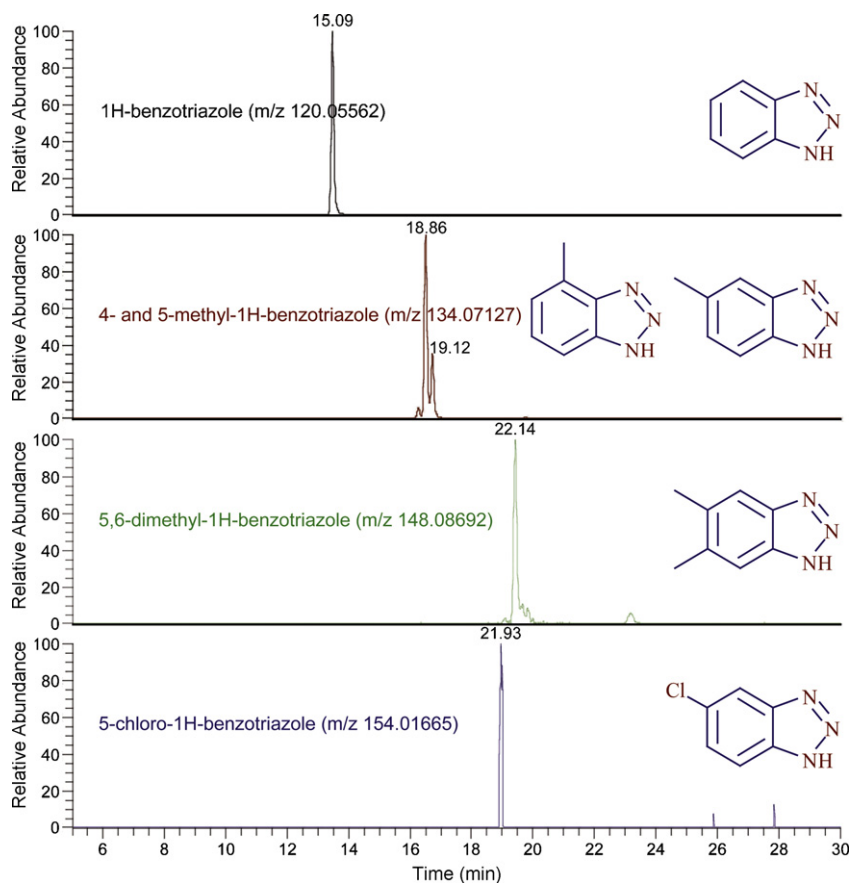


Fig. 1. LC-MS extracted-ion chromatograms (positive ion mode, mass window 7 ppm) of five benzotriazoles in a drinking water sample with concentrations ranging from <0.01 to 0.1  $\mu\text{g/L}$ .

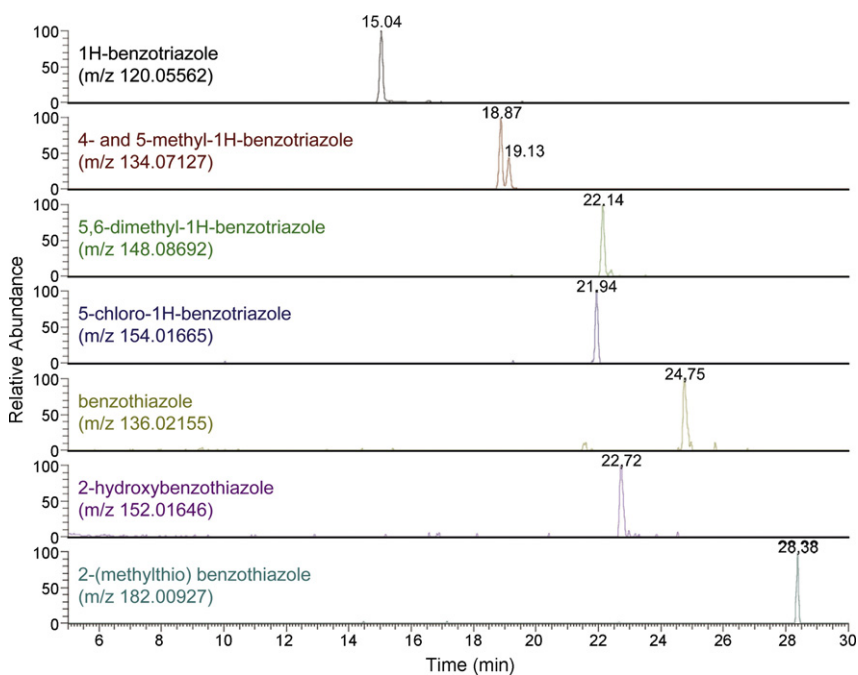


Fig. 2. LC-MS extracted-ion chromatograms (positive ion mode, mass window 7 ppm) of five benzotriazoles and three benzothiazoles in a surface water sample with concentrations ranging from <0.01 to 0.5  $\mu\text{g/L}$ .

In surface waters, eight out of ten compounds were detected (see Fig. 2), with concentrations ranging between 0.1 and 1.5  $\mu\text{g/L}$ . The compounds 1-hydroxybenzotriazole and 2-aminobenzothiazole were not detected in surface water.

In the effluents of both sewage treatment plants, eight out of ten compounds tested were detected and the maximum concentrations were observed for 1H-benzotriazole (8  $\mu\text{g/L}$ ) and methyl-1H-benzotriazole (summed concentration of two isomers: 3  $\mu\text{g/L}$ ).

#### 4. Conclusions

In Dutch drinking water samples, the compounds 1H-benzotriazole, 4- and 5-methyl-1H-benzotriazole, 5,6-dimethyl-1H-benzotriazole, 5-chloro-1H-benzotriazole and benzothiazole were detected. The concentration levels ranged from 0.01 to 0.2  $\mu\text{g/L}$ . In surface waters, eight out of ten compounds tested were actually observed in concentration levels ranging between 0.1 and 1.0  $\mu\text{g/L}$ . In addition, in the effluents of two sewage treatment plants, eight out of ten compounds tested were detected with maximum concentrations for 1H-benzotriazole of 8  $\mu\text{g/L}$  and for methyl-1H-benzotriazole of 3  $\mu\text{g/L}$  (summed concentration of two isomers). It can be concluded that 1H-benzotriazoles and benzothiazoles are omnipresent polar contaminants in the Dutch water cycle.

High resolution mass spectrometry coupled to LC is a very powerful combination for screening purposes. The application of

accurate mass screening described in this article demonstrates that current day analytical instrumentation is well equipped to meet the challenges posed by newly emerging polar chemicals.

#### References

- [1] S. Weiss, T. Reemtsma, *Anal. Chem.* 77 (2005) 7415.
- [2] A. Kloepper, M. Jekel, T. Reemtsma, *J. Chromatogr. A* 1058 (2004) 81.
- [3] A. Kloepper, M. Jekel, T. Reemtsma, *Environ. Sci. Technol.* 39 (2005) 3792.
- [4] D.A. Cancillia, J.C. Baird, R. Rosa, *Bull. Environ. Contam. Toxicol.* 70 (2003) 868.
- [5] W. Giger, C. Schaffner, D. Voutsas, *Proceedings of the 1st International Workshop on Liquid-Chromatography-Tandem Mass Spectrometry for Screening and Trace Level Quantitation in Environmental and Food Samples, Barcelona, Spain, 2005.*
- [6] S.D. Richardson, *Anal. Chem.* 78 (2006) 4021.
- [7] S.D. Richardson, *Anal. Chem.* 79 (2007) 4295.
- [8] A.C. Hogenboom, J.A. van Leerdam, P. de Voogt, *J. Chromatogr. A* 1216 (2009) 510.
- [9] B.K. Matuszewski, M.L. Constanzer, C.M. Chavez-Eng, *Anal. Chem.* 75 (2003) 3019.
- [10] Commission Decision, 2002/657/EC of 12 August 2002 implementing Council Directive 96/23/EC concerning performance of analytical methods and the interpretation of results, *Off. J. Eur. Commun.*, 2002, L221/8.
- [11] M.W.F. Nielen, M.C. van Engelen, R. Zuiderent, R. Ramaker, *Anal. Chim. Acta* 586 (2007) 122.
- [12] M. Krauss, J. Hollender, *Anal. Chem.* 80 (2008) 834.
- [13] D.A. Pillard, J.S. Cornell, D.L. Dufresne, M.T. Hernandez, *Water Res.* 35 (2001) 557.
- [14] ACD Laboratories, [www.acdlabs.com](http://www.acdlabs.com).